


The **ADONE** couch is a three-section couch. The height can be adjusted electrically. The couch is suitable for examination and traction therapy and is used by rehabilitation and physiotherapy professionals. The couch is equipped for lumbar and cervical traction and with a removable traction unit mounting and with retractable castors.



**TRACTION UNIT
NOT SUPPLIED**

TECHNICAL SPECIFICATIONS

Sections	3
Middle section	Fixed
Main section sliding by	12 cm
Under bed clearance	12 cm
Lateral offset	32 cm
Width	70 cm
Length	200 cm
Max working load, evenly distributed	200 kg
Height – electrically adjusted	55 - 86 cm
Power supply	230 V 50/60 Hz 1,5A
Actuator thrust	8000 N
Protection	IP 66
Headrest angle – with gas spring	+60° – 20°
Breathing hole and bung	as standard
Retractable castors	4
Weight	70 kg

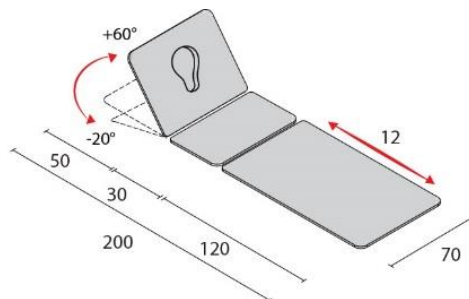

FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

ADONE COUCH MODELS

Description	Code
ADONE with foot pedal	114011

HEIGHT ADJUSTMENT SYSTEMS

ADONE with 24V foot pedal



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

TECHNICAL SPECIFICATIONS – PADDING

Name of material	Polyurethane foam
Reaction to fire class	1 IM (highest)
Density	30 kg/m ³
Height	40 mm



AVAILABLE UPHOLSTERY COLOURS

When placing your order, please specify the upholstery colour.



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

TECHNICAL SPECIFICATIONS – UPHOLSTERY

TECHNICAL DATA SHEET HAPPY AE I

<i>Resistance to penetration by body fluids (blood and urine)</i>	0-20 kPa	ISO 16603: 2004	NO PENETRATION
<i>Reaction ignition of upholstered furniture exposed to a small flame</i>		UNI 9175: 1987 UNI 9175/FA1: 1994	CLASS 1.IM <i>Item tested with flame retardant foam, density: 25-30 Kg/m³. Different foams can give different results.</i>
<i>Smouldering cigarette test -Propane flame test</i>		IMO - 2010 FTP CODE Part 8 "Test for upholstered furniture" par 7.2 - 7.3	PASS
<i>Specific features</i>	<i>Antibacterial</i>	<i>Fungicide</i>	
<i>Contains biocidal active substances Tiabendazolo and Folpet to help preservation of the surface and oppose degradation caused by microbial growth (according to EU Regulation n. 528/2012)</i>			
<i>Kindly note: pearl colours not suitable if used under intense crackling</i>			
<i>Composition</i>	<i>PLASTICISED P.V.C.</i>		88%
	<i>POLYESTER</i>		12%

Cleaning rules

Clean with damp soapy cloth and rinse well with clean water. Do not use solvents, bleaches, synthetic detergents, polishes or aerosol spray.

Attention

Normally synthetic leather light colours should not come in contact with no fixed closed clothes (such as jeans and similar) not to cause stains and rings that cannot be cleaned up.

The overmentioned data correspond to average values communicated in good faith, but they cannot be considered as a contract document. The product use has to be considered in terms of what has to be realized and it cannot be held up responsible when it is improperly used. Slight colour differences among several lots are admitted.

	STANDARDS	U.M.	MINIMUM ENVIRONMENTAL CRITERIA LIMITS	OBTAINED VALUES	
<i>Phthalates</i>	UNI EN ISO 14389: 2014	%	0,1	< 0,06	
<i>Arylamines</i>	UNI EN 14362-1: 2012	mg/Kg	≤ 30	NOT DETECTABLE	
<i>Determination of formaldehyde free and hydrolysed formaldehyde</i>	UNI EN ISO 14184-1:2011	mg/Kg	≤ 75	< 3	
<i>Determination of formaldehyde free and hydrolysed formaldehyde (classroom furniture)</i>	UNI EN ISO 14184-1:2011	mg/Kg	≤ 20	< 3	
<i>Determination of metal content</i>	UNI EN 16711-2:2015	mg/Kg	ANTIMONY	≤ 30	NOT DETECTABLE
			ARSENIC	≤ 1	NOT DETECTABLE
			CADMIUM	≤ 0,1	NOT DETECTABLE
			CHROME	≤ 2	NOT DETECTABLE
			COBALT	≤ 4	NOT DETECTABLE
			COPPER	≤ 50	NOT DETECTABLE
			LEAD	≤ 1	NOT DETECTABLE
			MERCURY	≤ 0,02	NOT DETECTABLE
			NICKEL	≤ 1	NOT DETECTABLE

* Test certificates available on request

FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

ADONE COUCH FULL EQUIPMENT

4 Ø 75 cm retractable castors for couch movement



Spreader bar for pelvic harness.



Thoracic harness.
Available in two sizes (small and large)



Pelvic harness.
Available in two sizes (small and large)



Head halter and spread bar



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

ADONE COUCH FULL EQUIPMENT

Trapezoid cushion. Dimensions: 40X50X30 cm.
The cushion colour corresponds with the couch colour.



Wedge cushion. Dimensions: 45x45x15 cm.
The cushion colour corresponds with the couch colour.



Half-roll cushion. Dimensions: 30x16x8 cm.
The cushion colour corresponds with the couch colour.



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956

GENERIC OPTIONAL ACCESSORIES

Code 120020 FLOOR PAPER ROLL HOLDER

Paper roll holder for examination, massage therapy, and treatment couches.

Easy to move, powder-coated metal structure.

Dimensions: length 71 cm, width 25 cm, height 81 cm



Code 120060 DOUBLE-PLY PAPER ROLL

6 rolls

Dimensions: width 60 cm – length 65 m.



FISIOTECH SRL

Operating Headquarters: Via Toniolo, 44/E-F - 31028 Vazzola (TV) - ITALY

Registered Office: Via Valdoni, 19 - 31013 Cimavilla di Codognè (TV) - ITALY

Tel. +39 0438470342 - Fax +39 0438478762 - mail : info@fisiotech.com - pec : fisiotech@emailsec.it

VAT Number IT 03765100262 - Rea TV 296956